



Received on 03 July 2018; received in revised form, 05 September 2018; accepted, 15 September 2018; published 01 March 2019

ANALYTICAL METHOD DEVELOPMENT AND VALIDATION FOR SIMULTANEOUS ESTIMATION OF HALOBETASOL PROPIONATE AND MUPIROCIN IN THE RATIO 1:40 BY U.V. SPECTROSCOPY AND RP- HPLC METHOD

Arpit Bana^{*}, Mansi A. Sathe and Sadhana J. Rajput

Quality Assurance Laboratory, Centre of Relevance and Excellence in Novel Drug Delivery Systems, Pharmacy Department, The Maharaja Sayajirao University of Baroda, Fatehgunj, Vadodara - 390002, Gujarat, India.

Keywords:

Halobetasol propionate,
Mupirocin, U.V Spectroscopy,
Anderson-Darling normality test,
RP-HPLC

Correspondence to Author:

Arpit Bana

Research Scholar,
Quality Assurance Laboratory,
Centre of Relevance and Excellence
in Novel Drug Delivery Systems,
Pharmacy Department, Shri G. H.
Patel Building, Donor's Plaza, The
Maharaja Sayajirao University of
Baroda, Fatehgunj, Vadodara -
390002, Gujarat, India.

E-mail: arpitbana94321@gmail.com

ABSTRACT: The present work encompasses the development of two classical UV spectroscopy methods and an RP-HPLC method for the concurrent assessment of Halobetasol propionate and Mupirocin inculcating a statistical approach. The UV spectroscopy methods developed herein include absorption ratio method and first derivative spectroscopy method. The wavelength maxima selected for Absorption ratio method were 240 and 220 nm for Halobetasol Propionate and Mupirocin respectively. The zero crossing points selected were 240.626 nm and 220.180 nm for Halobetasol Propionate and Mupirocin respectively for first derivative spectroscopy. The chromatographic separation was achieved by using CHROMBUDGET C18 (250 × 4.6 mm) 5 μm column, mobile phase consisting of Acetonitrile: phosphate buffer (65:35 v/v, pH 3.2), at 1 ml/min flow rate and detection wavelength was 230 nm. The retention time of Halobetasol propionate, and Mupirocin was found to be 8.647 ± 0.06 min and 3.357±0.123 min, respectively. Linearity was observed in the concentration range of 5.125 - 5.75 μg/ml for Halobetasol propionate and 5-30 μg/ml for Mupirocin respectively. The analytical method was validated as per ICH, Q2 (R1) guidelines. The developed method was simple, specific and economic, which can be used for simultaneous estimation of Halobetasol propionate and Mupirocin in the gel dosage form.

INTRODUCTION: Halobetasol Propionate:

Halobetasol propionate is a topical corticosteroid drug with molecular weight of 484.965g/mol and pKa value of 12.46 chemically it is (6S, 8S, 9R, 10S, 11S, 13S, 14S, 16S, 17R)-17-(2-chloroacetyl)-6,9-difluoro-11,17-dihydroxy-10, 13, 16-trimethyl-6, 7, 8, 11, 12, 14, 15, 16-octahydrocyclopenta[a]phenanthren-3-one **Fig. 1**.

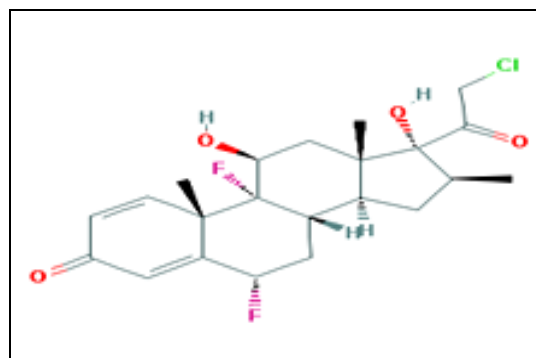


FIG. 1: STRUCTURE OF HALOBETASOL PROPIONATE

Halobetasol propionate is a potent corticosteroid that diffuses across cell membranes to interact with cytoplasmic corticosteroid receptors located in both the dermal and intra-dermal cells, thereby

QUICK RESPONSE CODE 	DOI: 10.13040/IJPSR.0975-8232.10(3).1392-01
	The article can be accessed online on www.ijpsr.com
DOI link: http://dx.doi.org/10.13040/IJPSR.0975-8232.10(3).1392-01	

activating gene expression of anti-inflammatory proteins mediated *via* the corticosteroid receptor response element. Halobetasol is mainly used for the treatment of edema, erythema, and pruritus through its cutaneous effect on vascular dilation and permeability¹.

Mupirocin: Mupirocin is a topical antibiotic with a molecular weight of 500.62g/mol and pKa value of 4.83. Chemically it is 9-[(E)-4-[(2S, 3R, 4R, 5S)-3,4-dihydroxy-5-[[[(2S, 3S)-3-[(2S, 3S) hydroxylbutan-2-yl] oxiran-2-yl] methyl] oxan-2-yl]-3-methylbut-2-enoyl]oxynonanoic acid **Fig. 2**. Mupirocin is a natural crotonic acid that is extracted from a strain of *Pseudomonas fluorescens*. Mupirocin inhibits bacterial protein synthesis by specific reversible binding to bacterial isoleucyl tRNA synthase with excellent activity against gram-positive staphylococci and streptococci; it is primarily used for the treatment of primary and secondary skin disorders, nasal infections, and wound healing (NCI04)¹.

No analytical methods are reported for this combination (Halobetasol propionate and mupirocin) in gel formulation however; analytical methods of single drug as well as with other combination are reported like, a spectrophotometric determination of clobetasol propionate, Halobetasol propionate, quinagolide hydrochloride, through charge transfer complexation², simultaneous determination of halobetasol propionate and fusidic acid-related substances by reversed phase high performance liquid chromatographic method³, Development and validation of RP-HPLC Method for simultaneous estimation Prednicarbate, Mupirocin and Ketoconazole in topical dosage form⁴, spectroscopic tools and implementation strategies for the chemical and pharmaceutical industries⁵, Formulation and evaluation of novel combined halobetasol propionate and fusidic acid ointment⁶, Practical High Performance Liquid Chromatographic method development⁷, Development and validation of a liquid chromatographic method for *in-vitro* mupirocin quantification in both skin layers⁸.

In this research article we have developed 2 UV methods for this combination A) Q-Absorbance ratio method B) First derivative zero crossing point method and an isocratic RP-HPLC method, we

confirm that this UV method and RP-HPLC method are not reported anywhere.

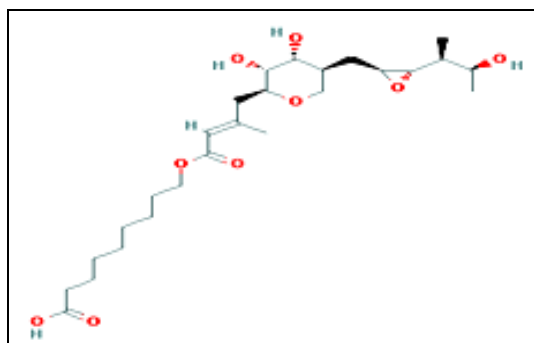


FIG. 2: STRUCTURE OF MUPIROCIN

MATERIALS AND METHODS:

Reagents and Chemicals: West-coast Pharma, Ahmedabad kindly gifted gift samples for Halobetasol propionate (HAL) and Mupirocin (MUP) API. The above-stated combination is an approved fixed-dose combination (FDC) by CDCSCO. The synthetic mixture of gel formulation was prepared in the laboratory on a pilot scale containing HAL 0.05% w/w and MUP 2.0% w/w **Table 1**.

Methanol analytical reagent grade (Fischer Scientific Pvt. Ltd., Mumbai, India) was used as the solvent for the UV spectrophotometric method. For chromatographic analysis, Acetonitrile of HPLC grade was procured from Rankem, India and used for mobile phase preparation. Double distilled water was prepared at the laboratory premises.

TABLE 1: COMPOSITION OF SYNTHETIC MIXTURE

Composition	Use	Strength in % w/w
Water	Excipient	48.55%
Disodium edetate	Stabilizer	0.05%
Ascorbic acid	Stabilizer	0.05%
Carbopol 934P	Thickening agent	1.0%
Poloxamer 407	Surfactant	0.2%
PEG-400	Co-Solvent	45.0%
Polysorbate 40	Surfactant	0.2%
Propylene glycol	Co-solvent	2.0%
BHA	Preservative	0.05%
BHT	Preservative	0.05%
Triethanolamine	Neutralizer	0.8%
Halobetasol Propionate	Drug	0.05%
Mupirocin	Drug	2.0%

Instrument and Software: Shimadzu UV-1700 double beam spectrophotometer connected to a computer loaded with Shimadzu UV probe 2.10 software was used for all the spectrophotometric measurements. Shimadzu UV-1800 double beam

spectrophotometer was also employed for ruggedness study. 1 cm quartz cells were used to measure the absorbance spectra of the reference and test solutions over the range of 200-400 nm. All the samples were weighed on an electronic analytical balance (A × 120, Shimadzu).

Preparation of Standard Stock Solution: Both UV spectrophotometric methods and RP-HPLC 10mg each of HAL and MUP were weighed accurately and transferred into a 10 ml volumetric flask to which methanol (Fischer Scientific Pvt. Ltd., Mumbai, India) was added up to the mark to produce a stock solution containing 1000 µg/ml of HAL and MUP respectively.

Selection of Analytical Wavelength: Solutions of 5 µg/ml of both drugs were prepared from working stock solution and scanned in the range of 200 nm to 400 nm against methanol as blank. The λ_{max} for Halobetasol propionate (HAL) and Mupirocin (MUP) was found to be 240 nm and 220 nm respectively and isoabsorptive point obtained was 230 nm for both drugs **Fig. 3**.

Method-A: Q-Absorbance Ratio Method:

- ✓ Absorbance ratio method uses the ratio of absorbance at two selected wavelengths, one which is an Iso-absorptive point and other being the λ_{max} of one of the two components.
- ✓ From the overlay spectra of two drugs, it is evident that HAL and MUP show an Iso-absorptive point at 230 nm **Fig. 3**. The second wavelength used is 220 nm, which is the λ_{max} of MUP.
- ✓ Six standard working solutions having concentration 5.125, 5.25, 5.375, 5.50, 5.625 and 5.75 µg/mL for HAL and 5, 10, 15, 20, 25, 30 µg/mL for MUP were prepared in methanol and the absorbance at 230 nm (isoabsorptive point) and 220 nm (λ_{max} of MUP) were measured and absorptivity coefficients were calculated.
- ✓ The absorbance of the sample solution (10 µg/ml of HAL and MUP), *i.e.* A_1 and A_2 were recorded at 230 nm (Iso-absorptive point) and 220 nm (λ_{max} of MUP) respectively, and ratios of absorbance were calculated, *i.e.*, A_2/A_1

- ✓ The relative concentration of two drugs in the sample was calculated using the following equations.

$$C_x = [(Q_m - Q_y) / (Q_x - Q_y)] \times A_1 / aX_1$$

$$C_y = [(Q_m - Q_x) / (Q_y - Q_x)] \times A_1 / aY_1$$

The Q-values and absorptivity for both drugs were calculated as follows,

$$Q_m = \frac{\text{Absorbance of Sample solution at 220 nm } (A_2)}{\text{Absorbance of Sample solution at 230 nm } (A_1)}$$

$$Q_x = \frac{\text{Absorptivity of MUP at 220 nm } (ax_2)}{\text{Absorptivity of MUP at 230nm } (ax_1)}$$

$$Q_y = \frac{\text{Absorptivity of HAL at 220nm } (ay_2)}{\text{Absorptivity of HAL at 230nm } (ay_1)}$$

Method-B: First Derivative Spectroscopy

Method: In this method, zero order spectra were derivatized to first order and then on basis of zero crossing points (ZCP) of the corresponding drugs the wavelength for analysis was chosen. As per the study, HAL was analyzed at 220.18 nm which was ZCP of MUP whereas MUP was analyzed at 240.626 nm which was ZCP of HAL. The parameters $\Delta\lambda$ is 5.0, and the scaling factor was set to be 20 **Fig. 4**.

Preparation of Calibration Curve of Halobetasol

Propionate: Concentration range for both drugs was decided according to the ratio (1:40) in which the combination of both drugs is available. For both UV spectrophotometric and RP-HPLC method, 0.25 ml (1000 µg/ml) of HAL standard was transferred into 25 ml volumetric flask to prepare 10µg/ml working standard of HAL from which aliquots of 0.125, 0.25, 0.375, 0.5, 0.625 and 0.75 ml were transferred into 10 ml volumetric flasks. A uniform spiking was carried out for all above samples in which 1.0 ml (50 µg/ml) of HAL standard was added to above volumetric flasks to obtain the final concentration range of 5.125-5.75 µg/ml.

Preparation of Calibration Curve of Mupirocin:

For MUP 2.5 ml (1000 µg/ml) of MUP was transferred into 25 ml volumetric flasks to prepare 100 µg/ml working standard from which aliquots of 0.5, 1.0, 1.5, 2.0, 2.5 and 3.0 ml were transferred into 10 ml volumetric flasks to obtain the final concentration range of 5-30 µg/ml of MUP.

Sample Preparation: Accurately-weighed gel equivalent to 50 mg of HAL and 2000 mg of MUP in 50 ml transferred to a volumetric flask, methanol was added and sonicated for 15 min to extract all the drugs from the gel formulation. The solution was cooled and then filtered by using 0.45 μ nylon filters. The filtrate was then evaporated and reconstituted with methanol and then transferred to 25 ml volumetric flask and diluted to the mark with methanol. The final concentration was achieved 5.5 μ g/ml for HAL and 20 μ g/ml for MUP. The mean values of peak areas of 5 such determinations were calculated, and the drug content in the gel was quantified using the calibration curve.

RP-HPLC Method Development:

Chromatographic Condition: Chromatographic separation was achieved by using Shimadzu LC-20AT high-performance liquid chromatography equipped with SPD-20A UV detector and CHROMBUDGET C18, (250 \times 4.6 mm ID, 5 μ m column) using Acetonitrile: Phosphate buffer (pH 3.2) (65:35 v/v) as mobile phase with flow rate 1.0 ml/min **Table 2**.

Preparation of Mobile Phase: Phosphate buffer, pH 3.2 was prepared by taking 0.2721 g of 20 mM potassium dihydrogen phosphate in a 100.0 ml volumetric flask, to which diluted orthophosphoric acid is added to bring the pH to acidic condition (pH 3.2). 35.0 ml of phosphate buffer (pH 3.2) and 65.0 ml of Acetonitrile were mixed and, sonicated for 15 min and filtered through 0.45 μ m membrane filters and used as mobile phase.

Selection of Detection Wavelength: Standard solutions of HAL and MUP were scanned between 200-400 nm in UV-visible spectrophotometer and showed good sensitivity at 230 nm as shown in **Fig. 5** which was selected as the detection wavelength.

Method Validation:

1. System Suitability Studies: System-suitability tests are an integral part of method development and are used to ensure adequate performance of the chromatographic system. Retention time (Rt), a number of theoretical plates (N) and tailing factor (T) were evaluated for six replicate injections of the drug at a concentration of 5 μ g/ml. The results which are given in **Table 3** were within acceptable limits. System suitability test is a pharmacopeia requirement.

2. Linearity: To evaluate the linearity, serial dilution of analyte were prepared from the stock solution was diluted with mobile phase to get a series of concentration ranging from 5-30 μ g/ml for MUP and 5.125-5.75 μ g/ml for HAL. From these solutions, 20 μ l injections of each concentration were injected into the HPLC system and chromatographed under the optimized conditions. A calibration curve was constructed by plotting the mean peak area (Y-axis) against the concentration (X-axis) and for UV methods A and B **Fig. 7-14** above concentration range of HAL and MUP were subjected to UV-Vis Spectrophotometer, and the calibration curve was plotted between Absorbance (Y-axis) and Concentration (X-axis)⁹ thus by this procedure regression and optical characteristics of UV and RP-HPLC method were obtained **Table 4** and **5**.

3. Limit of Detection (LOD): The detection limit of an individual analytical procedure is the lowest amount of analyte in a sample which can be detected but not necessarily quantified as an exact value. LOD was calculated using the following formula:

$$\text{LOD} = 3.3 * \sigma / S$$

Where σ is the standard deviation of the response and S is the slope of the calibration curve. The Limit of detection (LOD) for Halobetasol and Mupirocin is determined.

4. Limit of Quantification (LOQ): The quantification limit of an individual analytical procedure is the lowest amount of analyte in a sample which can be quantitatively determined with suitable precision and accuracy. LOQ was calculated using the following formula;

$$\text{LOQ} = 10 * \sigma / S$$

Where σ is the standard deviation of the response and S is the slope of the calibration curve. The limit of quantification for Halobetasol propionate and Mupirocin are determined⁹.

5. Precision: Intraday and Interday precision were measured in terms of % RSD. The experiment was repeated 3 times a day for intra-day and for 3 different days for inter-day precision. The average % RSD was found to be less than 2.0% for both UV methods **Table 6**.

To demonstrate agreement among results, a series of measurements were done with HAL and MUP. For RP-HPLC three replicate injections of the specific standard at various time intervals on the same day were injected into the system for intraday precision **Table 7** and were repeated on three different days for inter-day precision **Table 8**.

The % RSD (Relative Standard Deviation) of the results was calculated.

6. Accuracy: To check the Accuracy of different methods, Recovery studies were carried out from the pre-analyzed sample at three different levels of standard addition 80%, 100%, and 120%. Results of Recovery studies are shown in **Table 9** and **10** for both UV methods and **Table 11** for RP-HPLC method.

For each of the method explained above, % Recovery was the average of three determinations at each standard addition level. % Recovery for different methods was found to be between 98% - 102% which prove that all the methods were accurate.

7. Robustness: It should show the reliability of analysis concerning deliberate variations in method parameters. If measurements are susceptible to variation in analytical conditions, the analytical condition should be suitably controlled, or a precautionary statement should be included in the procedure. The evaluation of robustness should be considered during the development phase and depends on the type of procedure under study **Table 12**. It should show the reliability of an analysis concerning deliberate variations in method parameters. Robustness of the method was demonstrated by deliberately changing the chromatographic conditions like pH, mobile phase ratio and flow rate **Table 13**.

8. Ruggedness: Ruggedness means the ability of an analytical method to remain unaffected by small variations in method parameters and influential environmental factors and characterize its reliability during normal usage **Table 14** and **15**.

RESULT AND DISCUSSION:

Overlain UV Spectra of Halobetasol Propionate and Mupirocin for Q- Absorbance Ratio Method:

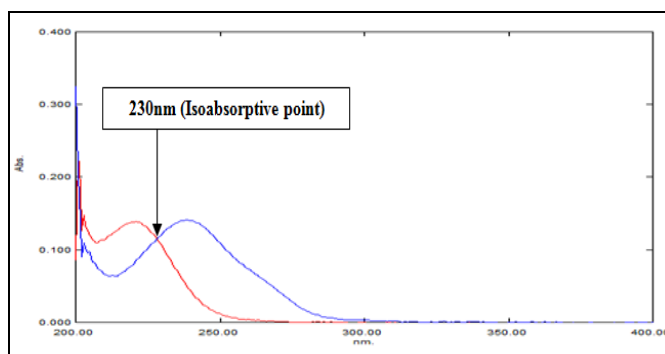


FIG. 3: OVERLAIN SPECTRA OF HAL AND MUP

Overlain UV Spectra of Halobetasol Propionate and Mupirocin for First Derivative ZCP Method:

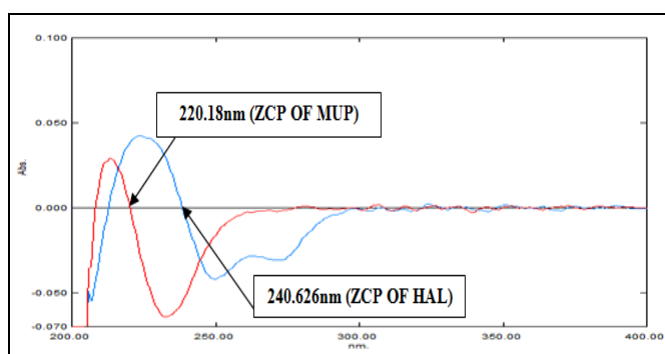


FIG. 4: OVERLAIN SPECTRA OF HAL (BLUE) AND MUP (RED) FOR FIRST DERIVATIVE ZCP METHOD

Selection of Detection Wavelength for RP-HPLC Method:

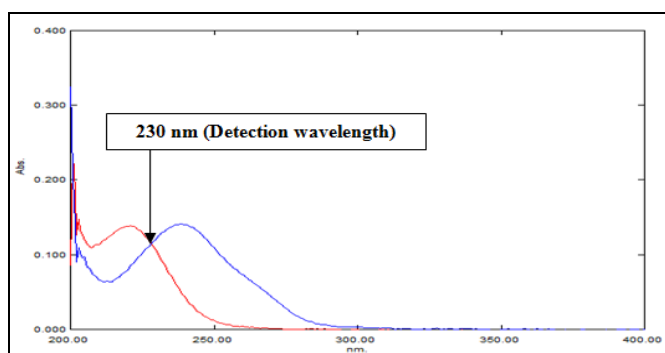


FIG. 5: OVERLAIN SPECTRA OF HALOBETASOL PROPIONATE AND MUPIROCIN

TABLE 2: OPTIMIZED CHROMATOGRAPHIC CONDITION

Method Parameter	Optimized Value
Column	CHROMBUDGET 100-5-C18 column, (Column dimensions: 250 mm × 4.6 mm, 5 μm)
Mobile Phase	ACN: Potassium Dihydrogen Phosphate, pH 3.2 = 65:35
Flow Rate	1.0 ml/min
Retention Time	3.357 min for MUP and 8.647 min for HAL
Detection Wavelength	230 nm
Temperature	Ambient

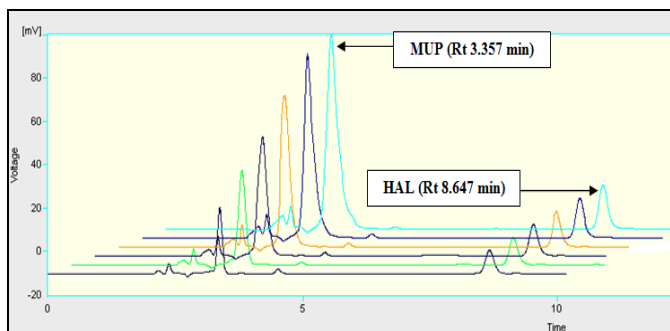


FIG. 6: 3D OVERLAIN CHROMATOGRAM OF HALOBETASOL PROPIONATE (RT 8.647 ± 0.06 min) AND MUPIROCIN (RT 3.357 ± 3.357 min)

TABLE 3: SST PARAMETERS

Parameter	Data Obtained	
	MUP	HAL
Retention time ± SD	3.357 ± 0.123	8.647 ± 0.06
Theoretical plate ± SD	4895 ± 18.98	12,323 ± 14.32
Tailing factor ± SD	0.895 ± 0.001	0.632 ± 0.043
Resolution ± SD	4.672 ± 0.501	18.181 ± 0.451

Calibration Curve of Halobetasol Propionate and Mupirocin for Q-Absorbance Ratio Method:

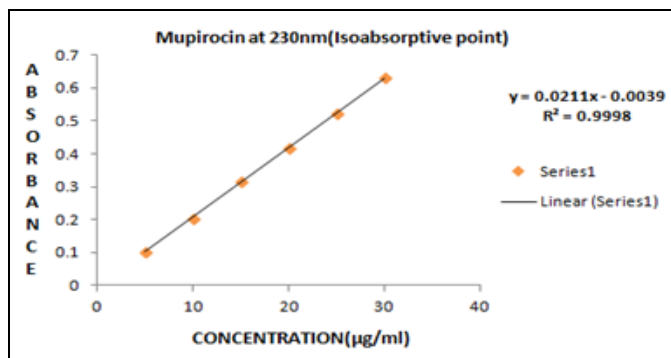


FIG. 7: CALIBRATION CURVE OF MUPIROCIN (MUP) AT 230 nm (ISOABSORPTIVE POINT)

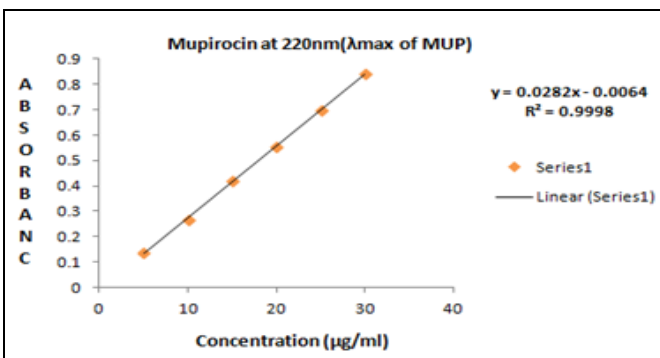


FIG. 8: CALIBRATION CURVE OF MUPIROCIN (MUP) AT 220 nm (λ_{MAX} OF MUP)

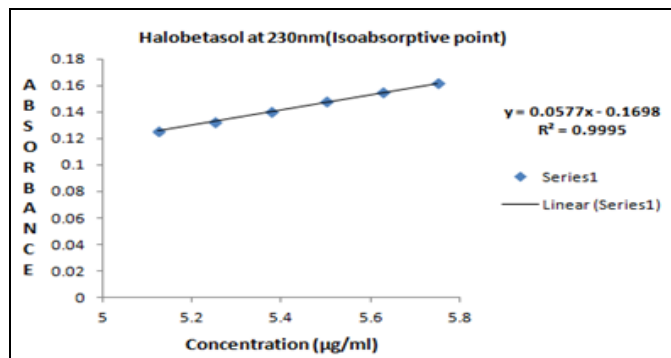


FIG. 9: CALIBRATION CURVE OF HALOBETASOL (HAL) AT 230 nm (ISOABSORPTIVE POINT)

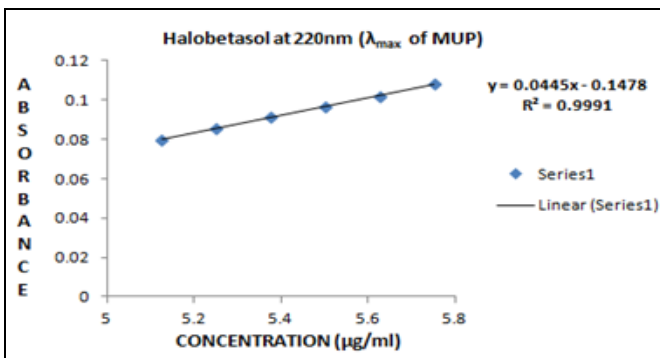


FIG. 10: CALIBRATION CURVE OF HALOBETASOL (HAL) AT 220 nm (λ_{MAX} OF MUP)

Calibration Curve of Halobetasol Propionate and Mupirocin for First Derivative ZCP Method:

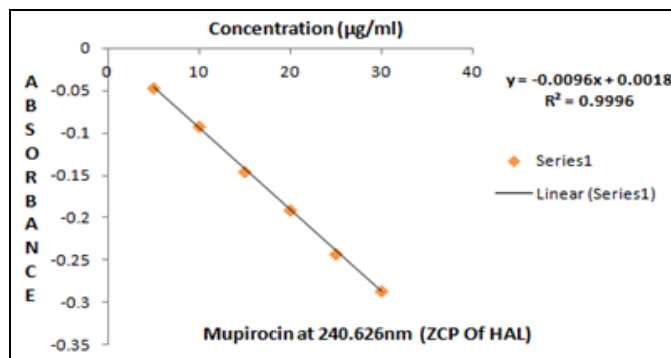


FIG. 11: CALIBRATION CURVE OF FIRST DERIVATIVE MUP AT 240.626 nm

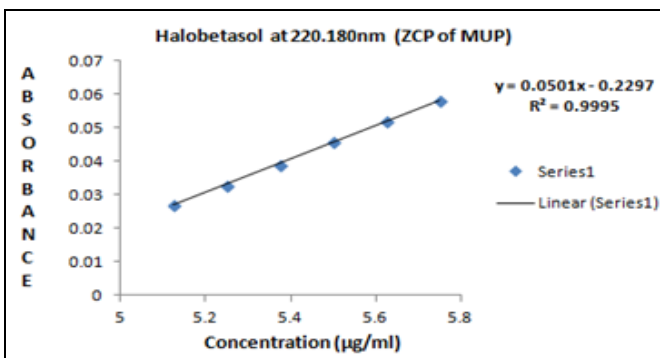
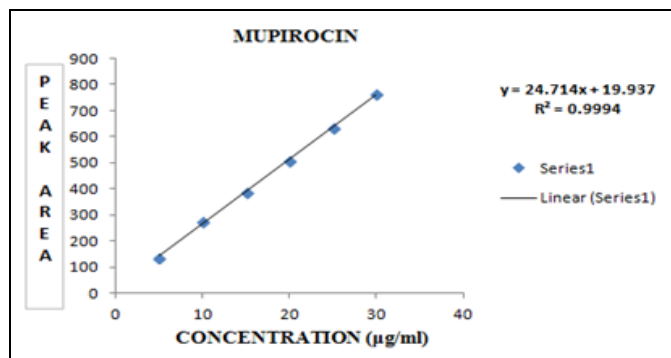
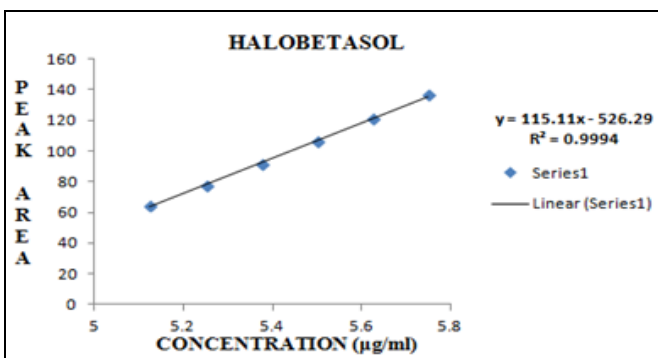


FIG. 12: CALIBRATION CURVE OF FIRST DERIVATIVE HAL AT 220.180 nm

Calibration Curve of Halobetasol Propionate and Mupirocin for Developed Isocratic RP-HPLC Method:**FIG. 13: CALIBRATION CURVE OF MUPIROCIN FOR RP-HPLC METHOD****FIG. 14: CALIBRATION CURVE OF HALOBETASOL FOR RP-HPLC METHOD****TABLE 4: REGRESSION AND OPTICAL CHARACTERISTICS FOR UV METHODS A & B**

Parameters	Q- Absorbance Ratio-method (A)				First Derivative ZCP method (B)	
	HAL		MUP		HAL	MUP
Analytical Wavelength (nm)	230 nm	220 nm	230 nm	220 nm	220.180 nm	240.626 nm
Beers's Range (µg/ml)	5.125-5.75		5-30		5.125-5.75	5-30
Slope	0.065	0.0445	0.0211	0.0282	0.0501	0.0096
Intercept	0.2098	0.1478	0.0039	0.0064	0.2297	0.0018
Correlation Coefficient	0.9994	0.9991	0.9998	0.9998	0.9995	0.9996
Standard Deviation of Intercept	0.075	0.0101	0.012093	0.000557	0.012093	0.000557
Limit of Detection (LOD) (µg/ml)	0.739	0.113	0.826	0.256	0.819	0.192
Limit of Quantification (LOQ) (µg/ml)	2.24	0.344	2.54	0.772	2.48	0.584

TABLE 5: REGRESSION AND OPTICAL CHARACTERISTICS OBTAINED FOR RP-HPLC METHOD

Parameter	Halobetasol	Mupirocin
Detection Wavelength (nm)	230 nm	
Beer's range (µg/ml)	5.125-5.75	5-30
Regression equation	Y=115.11X-526.29	Y=24.714X+19.93
LOD	0.2364	0.1069
LOQ	0.7095	0.3295

TABLE 6: INTRADAY AND INTERDAY PRECISION OF UV METHODS A & B

Parameter	Q-Absorbance Ratio-method (A)		First derivative ZCP method (B)	
	HAL	MUP	HAL	MUP
Intraday Precision (SD)	0.0013	0.0021	0.00044	0.00057
Intraday Precision (% RSD)	0.164	0.450	1.032	0.356
Interday Precision (SD)	0.00063	0.00105	0.00058	0.00103
Interday Precision (% RSD)	0.387	0.226	1.355	0.6494

TABLE 7: INTRADAY PRECISION OF RP-HPLC METHOD

Drug	Conc. (µg/ml)	Mean Area	Std. Dev	% RSD
Mupirocin	5	138.250	0.627	0.453
	15	390.900	3.893	0.995
	30	764.173	7.005	0.916
Halobetasol	5.125	66.527	0.530	0.796
	5.375	89.782	1.500	1.670
	5.75	140.249	1.792	1.277

TABLE 8: INTERDAY PRECISION OF RP-HPLC METHOD

Drug	Conc. (µg/ml)	Mean Area	Std. Dev	% RSD
Mupirocin	5	138.789	0.7514	0.541
	15	389.945	2.023	0.518
	30	765.145	7.058	0.922
Halobetasol	5.125	65.256	0.523	0.801
	5.375	90.065	1.022	1.134
	5.75	141.069	1.593	1.129

TABLE 9: RECOVERY DATA OF UV METHOD A

Recovery Drug	Q-Absorbance Ratio Method (A)			
	% Spiking	Conc. Added	Conc. Found	% Recovery \pm SD
Halobetasol	80	4.1	4.12	101.08% \pm 0.863
	100	5.125	5.14	101.21% \pm 1.12
	120	6.15	6.2	99.27% \pm 0.98
Mupirocin	80	4	4.25	101.04% \pm 0.98
	100	5	5.12	101.6% \pm 1.12
	120	6	5.96	99.44% \pm 0.892

TABLE 10: RECOVERY DATA OF UV METHOD B

Recovery Drug	First derivative ZCP method (B)			
	% Spiking	Conc. Added	Conc. Found	% Recovery \pm SD
Halobetasol	80	4.1	4.14	100.65% \pm 0.819
	100	5.125	5.128	100.86% \pm 0.921
	120	6.15	6.18	101.44% \pm 1.17
Mupirocin	80	4	3.95	99.00% \pm 0.89
	100	5	4.96	99.16% \pm 0.98
	120	6	5.96	99.55% \pm 0.93

TABLE 11: RECOVERY DATA OF RP-HPLC METHOD

% spiking	Concentration actual (μ g/ml)		Concentration added (μ g/ml)		Concentration recovered (μ g/ml)		% Recovery \pm standard deviation	
	HAL	MUP	HAL	MUP	HAL	MUP	HAL	MUP
80	5.125	5	4.1	4	4.14	3.95	100.6 \pm 0.22	100.41 \pm 0.76
100	5.125	5	5.125	5	5.128	4.96	100.25 \pm 0.232	101.3 \pm 0.22
120	5.125	5	6.15	6	6.18	5.96	99.79 \pm 0.85	99.8 \pm 0.01

TABLE 12: ROBUSTNESS DATA OF UV METHODS A & B

Brand of Methanol	Q-Absorbance Ratio Method (A)				First Derivative ZCP Method (B)			
	HAL	% RSD	MUP	% RSD	HAL	% RSD	MUP	% RSD
SPECTRO	0.144	0.68	0.136	0.491	0.027	0.83	0.051	1.17
	0.158		0.418		0.041		0.144	
	0.184		0.835		0.062		0.286	
RANKEM	0.146	1.69	0.134	0.452	0.028	1.75	0.049	1.05
	0.162		0.415		0.042		0.146	
	0.186		0.832		0.059		0.283	
FISCHER	0.143	0.66	0.135	0.584	0.028	1.66	0.049	0.536
	0.158		0.419		0.042		0.147	
	0.183		0.835		0.058		0.285	

TABLE 13: ROBUSTNESS DATA OF RP-HPLC METHOD

Factor	Retention Time (min)		Peak Area (mV.sec)	
	HAL	MUP	HAL	MUP
A. pH	3.1	3.350	8.510	139.026
	3.2	3.357	8.647	140.249
	3.3	3.410	8.695	141.116
	Mean \pm SD	3.372 \pm 0.025	8.617 \pm 0.0152	140.130 \pm 0.096
B. Flow Rate				
0.9ml/min	3.210	8.510	141.752	765.586
1ml/min	3.357	8.647	140.137	764.151
1.1ml/min	3.451	8.723	139.528	762.529
Mean \pm SD	3.33 \pm 0.050	8.626 \pm 0.0850	140.472 \pm 0.035	764.088 \pm 0.069
C. Buffer Ratio				
33	3.152	8.426	138.656	761.526
35	3.346	8.636	139.523	763.856
37	3.556	8.856	140.675	766.524
Mean \pm SD	3.351 \pm 0.052	8.639 \pm 0.015	139.618 \pm 0.01	763.968 \pm 0.108

TABLE 14: RUGGEDNESS DATA OF UV METHODS A & B

INSTRUMENT: UV-1800 SHIMADZU UV Spectrophotometer and UV-1700 SHIMADZU UV Spectrophotometer were used

Parameters	HAL		MUP	
	Method-A	Method-B	Method-A	Method-B
Mean (UV-1700&1800)	0.1617	0.042	0.4745	0.1598
Standard Deviation	0.00083	0.00049	0.00048	0.00070
%RSD	0.513	1.16	0.101	0.438

TABLE 15: RUGGEDNESS DATA OF RP-HPLC METHOD

ANALYST: Ruggedness of RP-HPLC Method was performed by using two different analyst

Parameters	HAL		MUP	
	Analyst-1	Analyst-2	Analyst-1	Analyst-2
Mean (Peak area)	140.130	140.472	763.968	764.324
Standard Deviation	0.096	0.112	0.108	0.128
%RSD	0.068	0.079	0.014	0.016

Anderson-Darling Normality Test: Anderson darling test for normality is one of three general normality tests designed to detect all departures from normality **Fig. 15** and **16**. The test rejects the hypothesis of normality when the p-value is less than or equal to 0.05. Failing the normality test allows to state with 95% confidence that given data

does not fit the normal distribution ¹⁰ **Table 16** and **17**.

$$AD = -N - \frac{2i-1}{N} (\ln(F(Y_i)) + \ln(1 - F(Y_{N+1-i})))$$

$$AD^* = AD \left(1 + \frac{0.75}{N} + \frac{2.25}{N^2} \right)$$

TABLE 16: ANDERSON DARLING DATA PLOT OF HALOBETASOL PROPIONATE

DATA	Sorted	Count	F1i	1-F1i	F2i	Si	N Plt Line	N Plt Line
0.05579	0.05579	1	0.08823154	0.91176846	0.0972061	-4.7587124	0.063695	0.1
0.11719	0.11719	2	0.20992205	0.79007795	0.1982854	-9.537201	0.135592	0.26
0.18127	0.18127	3	0.40599174	0.59400826	0.4009506	-9.076697	0.185323	0.42
0.23633	0.23633	4	0.59904936	0.40095064	0.5940083	-7.2329134	0.230811	0.58
0.30357	0.30357	5	0.80171463	0.19828537	0.7900779	-4.1096361	0.280542	0.74
0.35425	0.35425	6	0.90279388	0.09720612	0.9117685	-2.1409324	0.352438	0.9

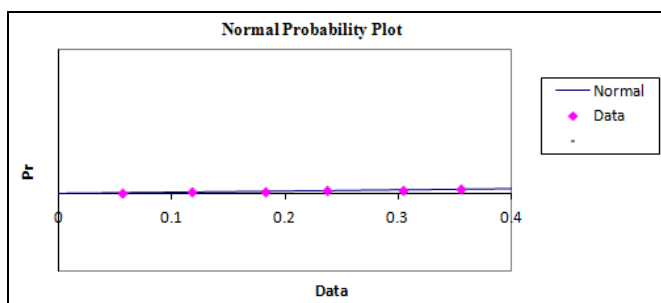


FIG. 15: NORMAL PROBABILITY PLOT OF HALOBETASOL PROPIONATE

6	Number of data points
0.208	Sample Mean
0.113	Sample Sigma
0.1427	AD test statistic
0.169435	AD* test statistic
0.934272	P-value

TABLE 17: ANDERSON DARLING PROBABILITY PLOT OF MUPIROCIN

DATA	Sorted	Count	F1i	1-F1i	F2i	Si	N Plt Line	N Plt Line
0.14111	0.14111	1	0.08855388	0.91144612	0.0946864	-4.7813291	0.160529	0.1
0.29041	0.29041	2	0.20449883	0.79550117	0.2102523	-9.4399199	0.34235	0.26
0.46301	0.46301	3	0.41300613	0.58699387	0.3908099	-9.1191337	0.468117	0.42
0.60461	0.60461	4	0.60919005	0.39080995	0.5869939	-7.1985612	0.583153	0.58
0.75513	0.75513	5	0.78974768	0.21025232	0.7955012	-4.1834227	0.70892	0.74
0.89954	0.89954	6	0.90531361	0.09468639	0.9114461	-2.1141633	0.890741	0.9

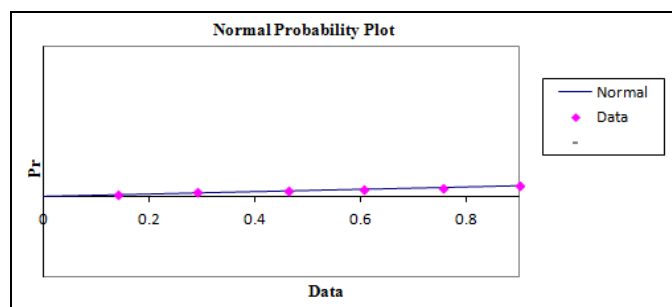


FIG. 16: NORMAL PROBABILITY PLOT OF MUPIROCIN

6	Number of data points
0.526	Sample Mean
0.285	Sample Sigma
0.1394	AD test statistic
0.165563	AD* test statistic
0.940611	P-value

CONCLUSION: Two UV methods Q-Absorbance ratio method & First derivative zero crossing point method and an RP- HPLC method were developed validated and applied for the simultaneous determination of Halobetasol Propionate and mupirocin in pharmaceutical gel formulation. The developed method was validated as per ICH guidelines and was found to be accurate, precise, robust and specific. The chromatographic elution step is undertaken in a short time. No, interference from any components of the pharmaceutical formulation was observed and the method has been successfully used to perform separation and identification of Halobetasol propionate and mupirocin.

ACKNOWLEDGEMENT: Authors like to thank West-coast pharmaceuticals, Ltd, Ahmedabad, Gujarat for providing free samples of Halobetasol propionate and mupirocin standards.

CONFLICT OF INTEREST: There is no conflict of interest between the authors and any other party

How to cite this article:

Bana A, Sathe MA and Rajput SJ: Analytical method development and validation for simultaneous estimation of Halobetasol propionate and Mupirocin in the ratio 1:40 by U.V. spectroscopy and RP- HPLC method. *Int J Pharm Sci & Res* 2019; 10(3): 1392-01. doi: 10.13040/IJPSR.0975-8232.10(3).1392-01.

and this article has not been published in any other journal.

REFERENCES:

1. Vanangamudi SS, Murali S and Madhavan S: Apex Labs Private Ltd, assignee. Methods and compositions for dermatological use comprising clobetasol and halobetasol and biopolymers. The United States patent application US 15/368,556. 2017.
2. Koo K, Jeon C and Bhutani T: Beyond monotherapy: a systematic review on creative strategies in topical therapy of psoriasis. *Journal of Dermatological Treatment* 2017; 28(8): 702-8.
3. Shrotriya S, Ranpise N, Satpute P and Vidhate B: Skin targeting of curcumin solid lipid nanoparticles-engrossed topical gel for the treatment of pigmentation and irritant contact dermatitis. *Artificial cells, nanomedicine, and biotechnology* 2018; 46(7): 1471-82.
4. Bassuoni YF, Elzanfaly ES, Essam HM and Zaazaa HE: Development and validation of stability indicating TLC densitometric and spectrophotometric methods for determination of clobetasol propionate. *Bulletin of Faculty of Pharmacy, Cairo University* 2016; 54(2): 165-74.
5. Bakeev KA: Process analytical technology: spectroscopic tools and implementation strategies for the chemical and pharmaceutical industries. John Wiley & Sons; 2010.
6. Curbete MM and Salgado HR: A critical review of the properties of fusidic acid and analytical methods for its determination. *Critical Reviews in Analytical Chemistry* 2016; 46(4): 352-60.
7. Praça FS, Medina WS, Eloy JO, Petrilli R, Campos PM, Ascenso A and Bentley MV: Evaluation of critical parameters for in vitro skin permeation and penetration studies using animal skin models. *European Journal of Pharmaceutical Sciences* 2018; 111: 121-32.
8. Cern A, Nativ-Roth E, Goldblum A and Barenholz Y: Effect of solubilizing agents on mupirocin loading into and release from PEGylated nanoliposomes. *Journal of pharmaceutical sciences* 2014; 103(7): 2131-8.
9. ICH I. Q2 (R1): Validation of analytical procedures: text and methodology. International Conference on Harmonization, Geneva 2005.
10. Rajput SJ and Sathe MA: Application of doe and statistical analysis for development and validation of analytical method for chlorhexidine gluconate and cetrimide in its bulk and pharmaceutical dosage forms. *International Journal of Pharmaceutical Sciences and Research* 2018; 9(7): 2800-6.